


Microrisk Lab

User Manual

Revised at January 2018

Version 1.0

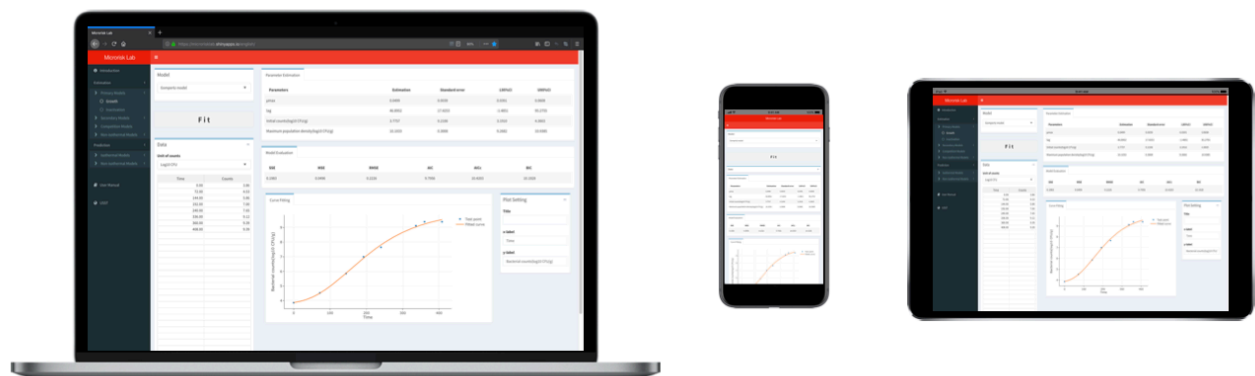

Disclaimer and Support

Disclaimer

Microrisk Lab and this manual provides NO WARRANTY. Microrisk Lab is copyrighted by University of Shanghai for Science and Technology (USST, Shanghai, China). It is not permitted to include Microrisk Lab in any other application. Besides, this tool is free to use but only for research purposes. We would very appreciate acknowledgement if the tool is used.

Feedback

If you have any suggestion or for technical questions for Microrisk Lab, please contact the developer and maintainer Yangtai Liu (usstlyt@163.com), Your comments are highly appreciated.

Document Revisions

| Date | Version | Document Changes |
|------------|----------|---------------------------|
| 01/09/2017 | Beta 0.1 | Initial Draft |
| 01/01/2018 | 1.0 | Draft for Updated version |

Contents

| | |
|---|-----------|
| Disclaimer and Support | 2 |
| 1 List of symbols | 4 |
| 2 Unit | 5 |
| 3 Programing basics | 5 |
| 4 Functions included in Microrisk Lab | 6 |
| 5 Layout of Microrisk Lab | 6 |
| 6 Estimation module of Microrisk Lab | 7 |
| <i>Practical example I - Isothermal growth fitting</i> | <i>7</i> |
| 7 Prediction module of Microrisk Lab | 10 |
| <i>Practical example II- Stochastic growth simulation</i> | <i>10</i> |
| 8 Predictive models integrated in Microrisk Lab | 14 |
| 9 Statistical indicators in Microrisk Lab | 17 |
| Reference | 18 |

1 List of symbols

| | |
|--------------------------------|--|
| $Y(t), Y_0, Y_{max}$ | the natural logarithm of real-time, initial, and maximum bacterial counts (ln CFU/g). |
| $y(t), y_0, y_{max}$ | the 10-base logarithm of real-time, initial, and maximum bacterial counts (log10 CFU/g). |
| y_{res} | the 10-base logarithm of the residual bacterial counts (log10 CFU/g). |
| μ_{max} | the maximum and optimal specific growth rate. |
| k_{max} | the maximum specific inactivation rate. |
| D | the time of decimal reduction. |
| D_{ref} | the reference time of decimal reduction. |
| t_{lag} | the time of lag. |
| S_l | the time of shoulder (or before inactivation). |
| t_{max} | the time arriving stationary phase in growth. |
| S_t | the time arriving tail (or stationary phase) in inactivation. |
| $T_{min}, T_{opt}, T_{max}$ | the minimum, optimal, and maximum temperature conditions for growth (°C). |
| T_{ref} | the reference temperature for inactivation(°C). |
| $pH_{min}, pH_{opt}, pH_{max}$ | the minimum, optimal, and maximum pH conditions for growth. |
| $aw_{min}, aw_{opt}, aw_{max}$ | the minimum, optimal, and maximum water activity conditions for growth. |
| q_0 | the initial physiological state of the bacterial in the Baranyi model. |
| δ, p | the coefficients in the Weibull model. |
| a, b | the coefficients in the square-root model. |
| A, m | the coefficients in the dynamic Huang model. |

2 Unit

The unit of bacterial count and time related variables can be defined by the user. The unit of predicted counting outputs will be transfer into 10-base logarithm. Note that the unit of the specific (growth/ inactivation) rate is a natural logarithm combined with unit of time, for example, ln CFU/g/h or ln CFU/g/min.

3 Programing basics

Microrisk Lab is developed by the open-source language R (version 3.5.1 for Mac OS X; <http://www.r-project.org>). All users are free to access and use this tool through the browser of any internet-connected device by the following link:

<http://microrisklab.shinyapps.io/english>

The operation of this Microrisk Lab must depend on the certain developed R packages, which were listed in Tab.1. All the required packages have been hosted and deployed in the Shinyapps.io sever (<https://www.shinyapps.io>).

Tab.1

| Package name | Version | Reference | Purpose |
|-----------------------|---------|-----------|--|
| <i>ggplot2</i> | 2.2.1 | | to generate visualized plots for output |
| <i>mc2d</i> | 0.1-18 | | to generate certain distribution for output |
| <i>Metrics</i> | 0.1.3 | | to calculate statistical indicators for output |
| <i>plotly</i> | 4.7.1 | | to generate interactive plots for output |
| <i>rhandsontable</i> | 0.3.6 | | to build interactive table for input |
| <i>shiny</i> | 1.0.5 | | to establish and upload the shiny app |
| <i>shinyalert</i> | 1.0 | | to pop the error alert for input and output |
| <i>shinydashboard</i> | 0.7.0 | | to build the interactive interface |
| <i>shinyWidgets</i> | 0.4.2 | | to build the interactive interface |
| <i>stats</i> | 3.4.3 | | to realize the regression analysis |

Microrisk Lab can be also used on computers without internet connection when installed locally. In this case, please contact the developer.

4 Functions included in Microrisk Lab

Microrisk Lab includes the following functions:

- Kinetic analysis of microbial isothermal growth
- Kinetic analysis of microbial non-isothermal growth
- Kinetic analysis of microbial isothermal inactivation
- Kinetic analysis of microbial non-isothermal inactivation
- Kinetic analysis of two-flora isothermal competition growth
- Secondary modeling of specific growth rate vs. temperature, pH and A_w .
- Deterministic/ Stochastic simulation for microbial isothermal growth
- Deterministic simulation for microbial isothermal growth
- Deterministic/ Stochastic simulation for microbial isothermal inactivation
- Deterministic simulation for microbial isothermal inactivation
- Output interactive plots of fitted and predicted curves.
- Output estimated results (estimates, standard error, and 95% confidential intervals) and multiple statistical indicators (SSE, MSE, RMSE, AIC, BIC) with respect to the experimental data in the 'Estimation' module.
- Output simulated bacterial counts or the distribution of the specific rate and final bacterial counts in the 'Prediction' module.
- Output correlation analysis between model parameters and simulated bacterial counts in the stochastic simulation.

5 Layout of Microrisk Lab

Fig.1 shows the page structure when loading in the Microrisk Lab via the browser in different devices. Users may switch the target question by the main menu in the left side. In the setting panel, user can input the experiment data and choose the model in here. The result panel will provide the estimated (predicted) values, statistical results, and interactive plots according to the setting.

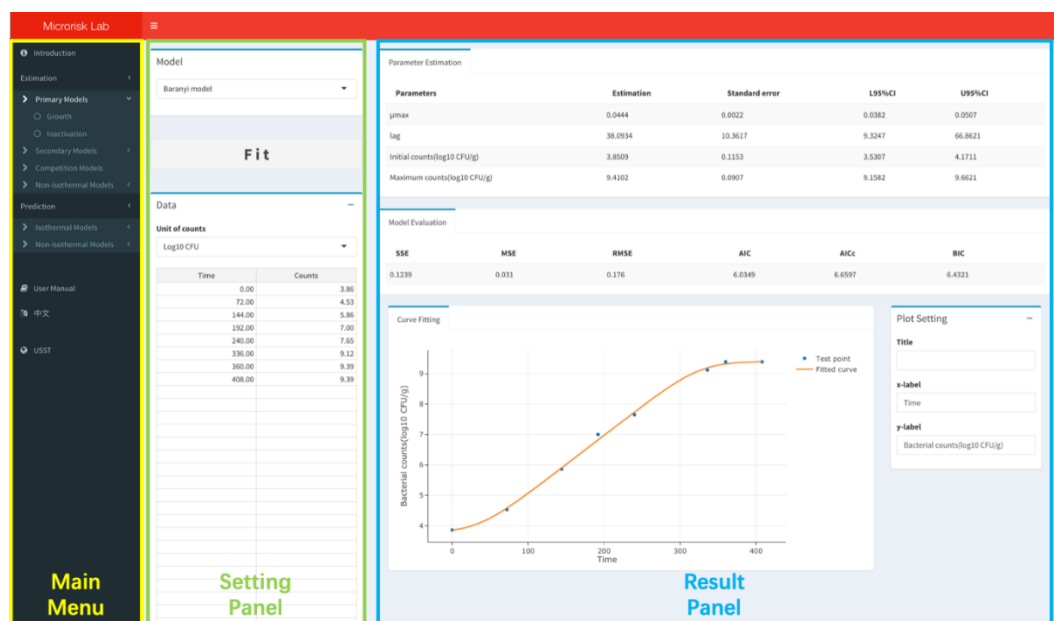


Fig.1 Typical layout of Microrisk Lab.

6 Estimation module of Microrisk Lab

The estimation module allows to solve multiple inverse problems in predictive microbiology, including ① isothermal growth fitting, ② isothermal inactivation fitting, ③④⑤ secondary model fitting, ⑥ two flora competition growth fitting, ⑦ non-isothermal growth fitting, and ⑧ non-isothermal inactivation fitting (Fig.2).

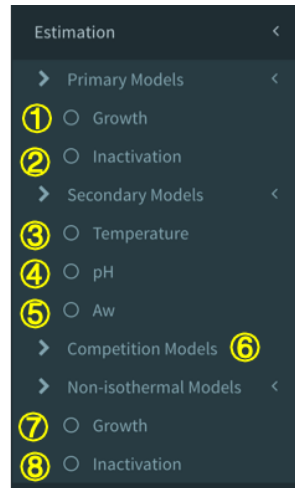


Fig.2 Different sections of model fitting in the estimation module.

Practical example I - Isothermal growth fitting

- (1) Choose ① the 'Growth' in the section of the 'Primary Models', and the setting panel of isothermal growth model will show up (Fig.3).

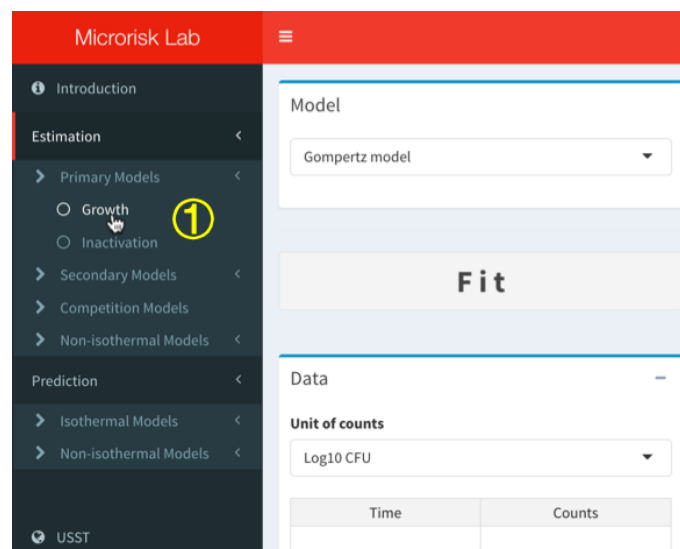


Fig.3 Layout of a section of the 'Estimation' module.

- (2) The experimental data can be ① directly typed (or ② copied from other table files) in the 'Data' box. Specifically, ③ the unit of bacterial counts should be confirmed by the user. If the inputted observations are more than 30, please ④ right click the mouse or ⑤ drag the last column to add additional columns. Here, a group of *Listera monocytogenes/ innocua* growth in tryptose phosphate broth (TPB) obtained from the ComBase database (www.combase.cc, ComBase ID: LM127_11) was used as the test dataset (Fig.4).

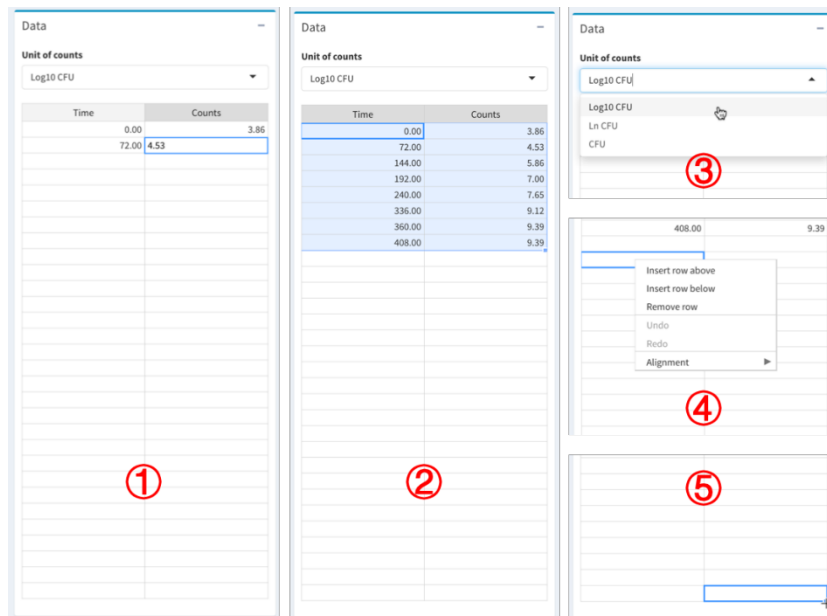


Fig.4 Boxes for the data input and unit selection.

- (3) After entering the data for model fitting, the growth model can be selected in the 'Model' list (Fig.5).

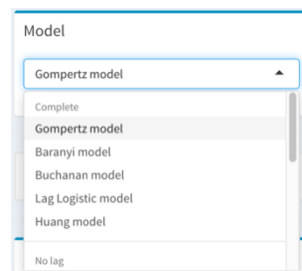


Fig.5 Box for model selection

- (4) Click ① the 'Fitting' button. After a necessary loading time, if the regression can be solved successfully, the ② estimated result, ③ evaluated result, and ④ interactive plot of the observation and fitted curve will show in the result panel (Fig.6A). Otherwise, a popup message will appear for the regression warning, which means that the non-linear regression is failed (Fig.6B). In this case, please check the unit of bacterial counts and the model selection.

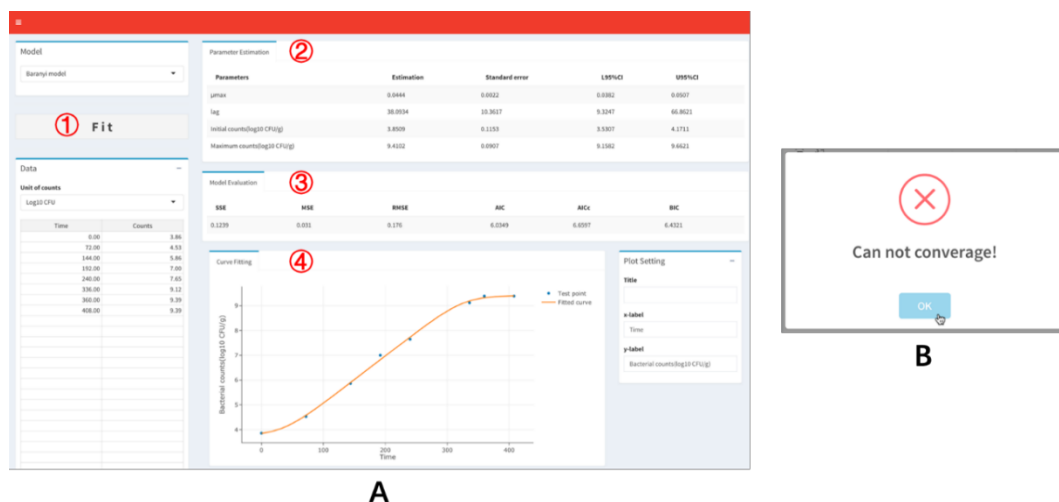


Fig.6 Layout of the interface after model fitting.

- (5) The observed and predicted value can be viewed on ① the interactive plot. The observed data or fitting curve can be omitted from the plot by clicking ② the legend. Meanwhile, it is easy to edit the axis detail (③ the range and ④ title) of the interactive plot in real-time by ⑤ the box of 'Plot Setting'. After all, the plot is adjustable and downloadable by using ⑥ the 'Plotly toolbox' (Fig.7).

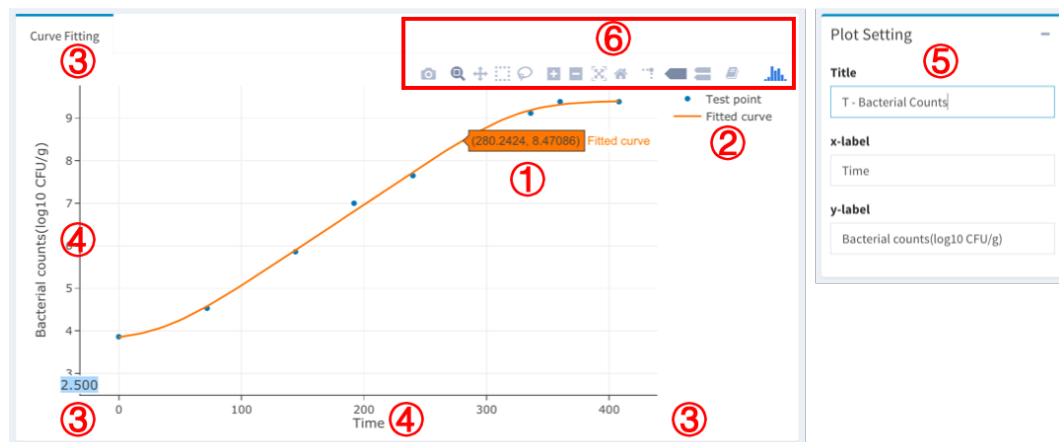


Fig.7 The interactive plot and the editorial box.

- (6) Note that, in the section of 'Non-isothermal Models', the additional time-temperature profile is needed to be imported for dynamic fitting.

7 Prediction module of Microrisk Lab

The prediction module allows to solve the ①② isothermal and ③④ non-isothermal forward problem in predictive microbiology (Fig.8). There are no limitations in the condition setting. Users may simulate the bacterial growth or inactivation with the prior knowledge on the kinetic parameter and growth/ death boundary. Moreover, both deterministic and stochastic models are provided in the isothermal simulation.

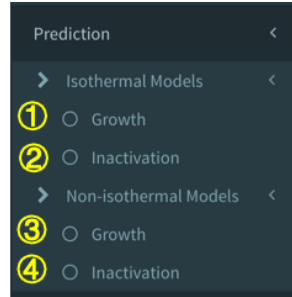


Fig.8 Different sections of model fitting in the estimation module.

Practical example II- Stochastic growth simulation

- (1) The condition setting of the growth simulation is adopted from the stochastic growth of *Salmonella* Typhimurium individual cells researched by Koutsoumanis and Lianou (2013). Tab.2 lists the setting for simulation. The Buchanan model is chosen as the growth model for individual cells. A 10,000 times iteration was realized based on the Monte-Carlo simulation method.

| Tab.2 Stochastic growth simulation settings for Microrisk Lab | | |
|---|--------------------|-----------|
| Parameters | Microrisk Lab | |
| y_0 (log ₁₀ CFU/g) | Distribution | Normal |
| | Mean | 0 |
| | Standard deviation | 0 |
| y_{max} (log ₁₀ CFU/g) | Distribution | Normal |
| | Mean | 8 |
| | Standard deviation | 0 |
| t_{lag} | Distribution | LogNormal |
| | Mean | 3.355 |
| | Standard deviation | 0.896 |
| | Shift | -1.628 |
| μ_{max} | Distribution | Logistic |
| | Mean | 0.754 |
| | Standard deviation | 0.024 |
| t | Distribution | Uniform |
| | Maximum | 0 |
| | Minimum | 8 |
| Model | Buchanan model | |
| Iteration times | 10,000 | |

- (2) Choose ① the 'Growth' section of the 'Isothermal Models' in the 'Prediction' module, and ② choose 'Stochastic' model type in the setting panel (Fig.9). Then set the ③ 'Iteration time' and ④ 'Model' to '10,000' and 'Buchanan model', respectively.

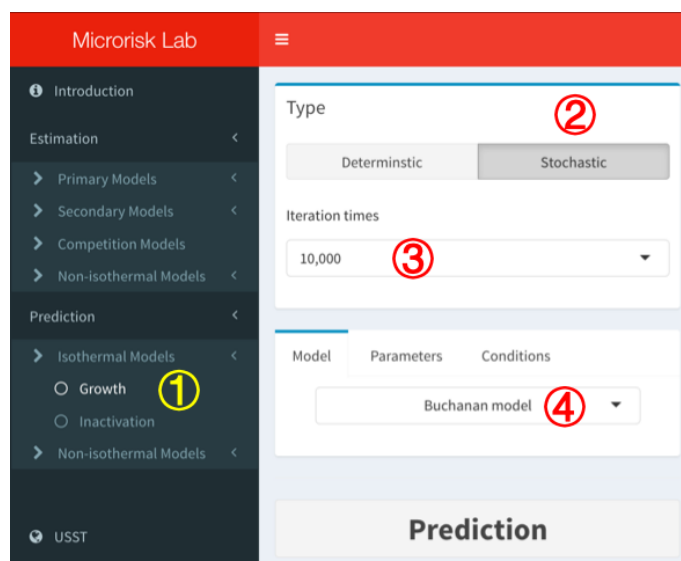


Fig.9 Layout of a section of the 'Prediction' module.

- (3) Switch to the ① 'Parameters' tab to determine the setting of the (distribution of) ② y_{max} , ③ t_{lag} , and ④ μ_{max} according to Tab.2 (Fig.10).

Fig.10 Box for kinetic parameter setting.

- (4) Switch to the ① ‘Conditions’ tab to determine the setting of the (distribution of) ② y_0 , ③ t according to Tab.2 (Fig.11).


Fig.11 Box for condition setting.

- (5) Click the ① ‘Prediction’ button. After a necessary loading time, if no contradiction in the setting, the ② simulated curve/point and ③ predicted result will show in the result panel (Fig.12A). Otherwise, different popup messages will appear for the simulation warning (Fig.12B-D). In these cases, please check the setting of kinetic parameters and the condition of simulated environment.




Fig.12 Layout of the interface after simulation.

- (6) The stochastic growth simulation can be viewed on ① the interactive plot, which is also adjustable and downloadable (Fig.13A). The distribution of ② the estimated μ_{max} and ③ final bacterial concentration (y_{final}), as well as ④ the estimated mean value and standard deviation will be presented and listed (Fig.13B). The sensitivity analysis on model parameters is realized by calculating the Pearson correlation between different factors and the bacterial counts. Here, according to ⑤ the correlation plot, the duration of growth time is the most sensitive parameter for the bacterial counts during the stochastic growth of *S. Typhimurium* single cell (Fig.13B).

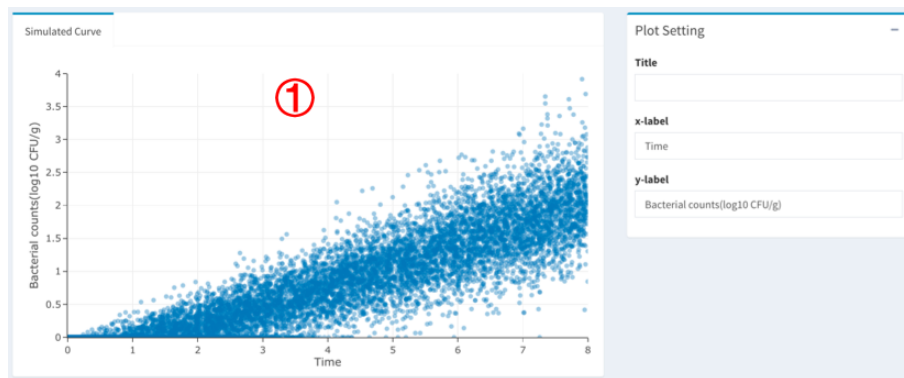

A

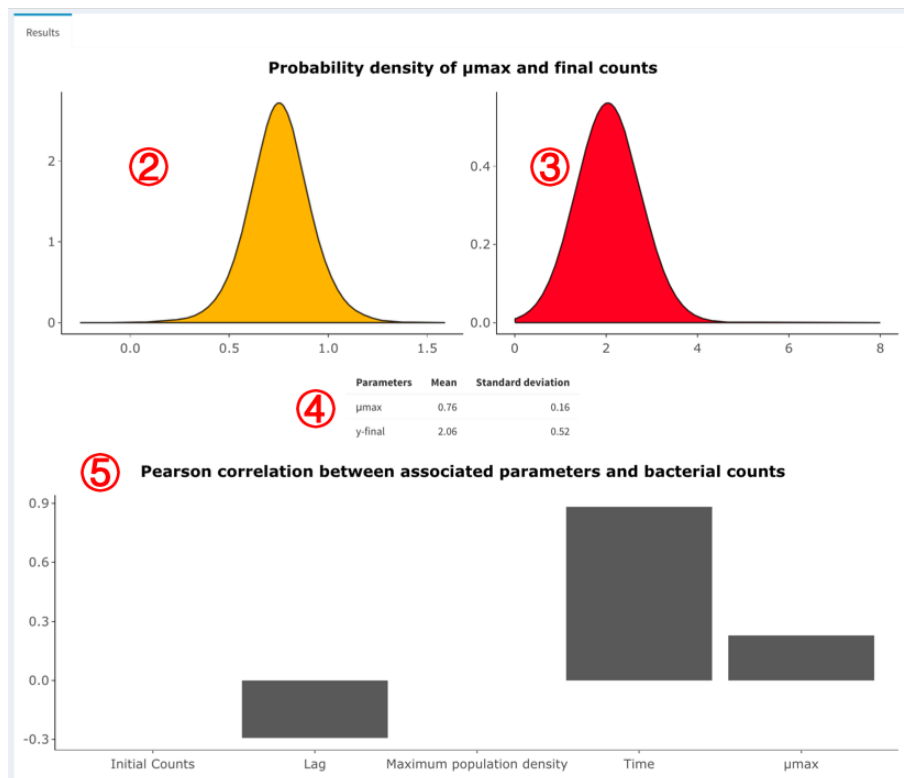

B

Fig.13 The result of stochastic simulation.

(7) Note that, in the section of 'Non-isothermal Models', only deterministic model is provided for now.

8 Predictive models integrated in Microrisk Lab

Microrisk Lab consists of 11 isothermal growth models (Tab.2), 9 inactivation models (Tab.3), 10 secondary models (Tab.4), 2 competition growth models (Tab.5), and 4 non-isothermal models (Tab.6) for estimation or prediction works.

Tab.2. Isothermal growth models included in Microrisk Lab

| Name | Formula |
|---------------------------------|--|
| Complete model | |
| Gompertz model ¹ | $Y(t) = Y_0 + (Y_{max} - Y_0) \exp \left\{ -\exp \left[\frac{2.71 \mu_{max} (t_{lag} - t)}{Y_{max} - Y_0} + 1 \right] \right\}$ |
| Baranyi model ² | $\begin{cases} Y(t) = Y_0 + A(t) - \ln \left[1 + \frac{\exp(\mu_{max} A(t)) - 1}{\exp(Y_{max} - Y_0)} \right] \\ A(t) = \mu_{max} t + \ln[\exp(-\mu_{max} t) + \exp(-h_0) - \exp(-\mu_{max} t - h_0)] \\ h_0 = \mu_{max} t_{lag} \end{cases}$ |
| Buchanan model ³ | $\begin{cases} y(t) = y_0, & t < t_{lag} \\ y(t) = y_0 + \frac{\mu_{max}}{\ln 10} (t - t_{lag}), & t_{lag} \leq t < t_{max} \\ y(t) = y_{max}, & t \geq t_{max} \end{cases}$ |
| Lag logistic model ⁴ | $\begin{cases} Y(t) = Y_0, & t < t_{lag} \\ Y(t) = Y_{max} - \ln \{ 1 + [\exp(Y_{max} - Y_0) - 1] \exp[-\mu_{max} (t - t_{lag})] \}, & t \geq t_{lag} \end{cases}$ |
| Huang model ⁵ | $\begin{cases} Y(t) = Y_0 + Y_{max} - \ln \{ \exp(Y_0) + [\exp(Y_{max}) - \exp(Y_0)] \exp(-\mu_{max} B(t)) \} \\ B(t) = t + \frac{1}{\alpha} \ln \frac{1 + \exp[-\alpha(t - t_{lag})]}{1 - \exp(\alpha t_{lag})} \end{cases}$ |
| No lag model | |
| Logistic model ⁶ | $Y(t) = Y_0 + Y_{max} - \ln \{ \exp(Y_0) + [\exp(Y_{max}) - \exp(Y_0)] \exp(-\mu_{max} t) \}$ |
| Buchanan model ⁷ | $\begin{cases} y(t) = y_0 + \frac{\mu_{max}}{\ln 10} t, & t < t_{max} \\ y(t) = y_{max}, & t \geq t_{max} \end{cases}$ |
| Reduced model | |
| Baranyi model ⁸ | $\begin{cases} Y(t) = Y_0 + \mu_{max} t + \ln[\exp(-\mu_{max} t) + \exp(-h_0) - \exp(-\mu_{max} t - h_0)] \\ h_0 = \mu_{max} t_{lag} \end{cases}$ |
| Buchanan model ⁹ | $\begin{cases} y(t) = y_0, & t < t_{lag} \\ y(t) = y_0 + \frac{\mu_{max}}{\ln 10} (t - t_{lag}), & t \geq t_{lag} \end{cases}$ |
| Huang model ¹⁰ | $Y(t) = Y_0 + \mu_{max} t + \frac{1}{\alpha} \mu_{max} \ln \frac{1 + \exp[-\alpha(t - t_{lag})]}{1 - \exp(\alpha t_{lag})}$ |
| Linear model | $Y(t) = Y_0 + \mu_{max} t$ |

¹Zwietering et al., 1990; ^{2/8} Baranyi and Roberts, 1995; ^{3/7/9} Buchanan et al., 1997; ⁴ Rosso et al., 1996; ^{5/6/10} Huang, 2008.

Tab.3. Isothermal inactivation models included in Microrisk Lab

| Name | Formula |
|--|--|
| Completed Geeraerd model ¹ | $y(t) = y_{res} + \log_{10} \left[\frac{(10^{y_0 - y_{res}} - 1) \exp(k_{max} S_l)}{\exp(k_{max} t) + \exp(k_{max} S_l) - 1} + 1 \right]$ |
| Three-phase model ² | $\begin{cases} y(t) = y_0, & t < S_l \\ y(t) = y_0 + \frac{k_{max}}{\ln 10} (t - S_l), & S_l \leq t < S_t \\ y(t) = y_{res}, & t \geq S_t \end{cases}$ |
| Weibull-tail model ³ | $y(t) = y_{res} + \log_{10} \left[(10^{y_0 - y_{res}} - 1) 10^{-\left(\frac{t}{\delta}\right)^p} + 1 \right]$ |
| No shoulder Geeraerd model ⁴ | $y(t) = y_{res} + \log_{10} \{ (10^{y_0 - y_{res}} - 1) \exp(k_{max} t) + 1 \}$ |
| No shoulder two-phase model ⁵ | $\begin{cases} y(t) = y_0 + \frac{k_{max}}{\ln 10} t, & t < S_t \\ y(t) = y_{res}, & t \geq S_t \end{cases}$ |
| No tail Geeraerd model ⁶ | $y(t) = y_0 + \frac{k_{max} t}{\ln 10} + \log_{10} \left\{ \frac{\exp(k_{max} S_l)}{1 + [\exp(k_{max} S_l) - 1] \exp(k_{max} t)} \right\}$ |
| No tail two-phase model ⁷ | $\begin{cases} y(t) = y_0, & t < S_l \\ y(t) = y_0 + \frac{k_{max}}{\ln 10} (t - S_l), & t \geq S_l \end{cases}$ |
| Weibull model ⁸ | $y(t) = y_0 - \left(\frac{t}{\delta}\right)^p$ |
| Bigelow model ⁹ | $y(t) = y_0 - \frac{t}{D}$ |

^{1/4/6} Geeraerd et al., 2000; ^{2/5/7} Buchanan and Golden, 1995; ³ Albert and Mafart, 2005; ⁸ van Boekel, 2002; ⁹ Bigelow, 1921.

Tab.4. Secondary models for μ_{max} included in Microrisk Lab

| Name | Formula |
|---|--|
| Temperature models | |
| Suboptimal square-root model ¹ | $\mu_{max} = [a(T - T_{min})]^2$ |
| Full square-root model ² | $\mu_{max} = \langle a(T - T_{min})\{1 - \exp[b(T - T_{max})]\} \rangle^2$ |
| Suboptimal Huang square-root model ³ | $\mu_{max} = [a(T - T_{min})^{0.75}]^2$ |
| Full Huang square-root model ⁴ | $\mu_{max} = \langle a(T - T_{min})^{0.75}\{1 - \exp[b(T - T_{max})]\} \rangle^2$ |
| Cardinal parameter model ⁵ | $\mu_{max} = \frac{\mu_{opt}(T - T_{max})(T - T_{min})^2}{[(T_{opt} - T_{min})(T - T_{opt}) - (T_{opt} - T_{max})(T_{opt} + T_{min} - 2T)](T_{opt} - T_{min})}$ |
| pH models | |
| Cardinal 3-parameter model ⁶ | $\mu_{max} = \frac{\mu_{opt}(pH - pH_{min})[pH - (2pH_{opt} - pH_{min})]}{(pH - pH_{min})[pH - (2pH_{opt} - pH_{min})] - (pH - pH_{opt})^2}$ |
| Cardinal 4-parameter model ⁷ | $\mu_{max} = \frac{\mu_{opt}(pH - pH_{min})(pH - pH_{max})}{(pH - pH_{min})(pH - pH_{max}) - (pH - pH_{opt})^2}$ |
| Quasi-mechanistic model ⁸ | $\mu_{max} = \mu_{opt}(1 - 10^{pH_{min} - pH})$ |
| Water activity models | |
| Cardinal 2-parameter model ⁹ | $\mu_{max} = \frac{\mu_{opt}(aw - aw_{min})^2}{(1 - aw_{min})^2}$ |
| Cardinal 3-parameter model ¹⁰ | $\mu_{max} = \frac{\mu_{opt}(aw - 1)(aw - aw_{min})^2}{(aw_{opt} - aw_{min})[(aw_{opt} - aw_{min})(aw - aw_{opt}) - (aw_{opt} - 1)(aw_{opt} + aw_{min} - 2aw)]}$ |

^{1/2} Ratkowsky et al., 1983; ^{3/4} Huang et al, 2011; ⁵ Rosso et al, 1993; ^{6/7} Rosso et al, 1995; ⁸ Presser et al. 1997; ^{9/10} Rosso and Robinson, 2001

Tab.5. Two flora competition growth models included in Microrisk Lab

| Name | Formula |
|--|--|
| Jameson - No lag Buchanan model ¹ | $\begin{cases} y_1(t) = \begin{cases} y_1 + \frac{\mu_{max1}}{\ln 10} t, & t < t_{max} \\ y_1 + \frac{\mu_{max1}}{\ln 10} t_{max}, & t \geq t_{max} \end{cases} \\ y_2(t) = \begin{cases} y_2 + \frac{\mu_{max2}}{\ln 10} t, & t < t_{max} \\ y_2 + \frac{\mu_{max2}}{\ln 10} t_{max}, & t \geq t_{max} \end{cases} \end{cases}$ |
| Jameson - Buchanan model ² | $\begin{cases} y_1(t) = \begin{cases} y_1, & t < t_{lag1} \\ y_1 + \frac{\mu_{max1}}{\ln 10} (t - t_{lag1}), & t_{lag1} \leq t < t_{max} \\ y_1 + \frac{\mu_{max1}}{\ln 10} (t_{max} - t_{lag1}), & t \geq t_{max} \end{cases} \\ y_2(t) = \begin{cases} y_2, & t < t_{lag2} \\ y_2 + \frac{\mu_{max2}}{\ln 10} (t - t_{lag2}), & t_{lag2} \leq t < t_{max} \\ y_2 + \frac{\mu_{max2}}{\ln 10} (t_{max} - t_{lag2}), & t \geq t_{max} \end{cases} \end{cases}$ |

* The inferior number 1 or 2 in competition growth models represent the flora type;

^{1/2} Vimont et al., 2006

Tab.6. Non-isothermal models included in Microrisk Lab

| Name | Formula |
|---|--|
| Non-isothermal growth models | |
| Baranyi - Cardinal parameter model ¹ | $\left\{ \begin{array}{l} \frac{dY(t)}{dt} = \mu_{max} \left[\frac{1}{1 + \exp(-Q(t))} \right] [1 - \exp(Y(t) - Y_{max})] \\ \frac{dQ(t)}{dt} = \mu_{max} \\ Q(t) = \ln \frac{q(t)}{1 - q(t)} \\ \mu_{max} = \frac{\mu_{opt}(T - T_{max})(T - T_{min})^2}{[(T_{opt} - T_{min})(T - T_{opt}) - (T_{opt} - T_{max})(T_{opt} + T_{min} - 2T)](T_{opt} - T_{min})} \end{array} \right.$ |
| Huang - Cardinal parameter model ^{2/3} | $\left\{ \begin{array}{l} \frac{dY(t)}{dt} = \mu_{max} \left[\frac{1}{1 + \exp(-4(t - t_{lag}))} \right] [1 - \exp(Y(t) - Y_{max})] \\ t_{lag} = \frac{\exp(A)}{\mu_{max}^m} \\ \mu_{max} = \frac{\mu_{opt}(T - T_{max})(T - T_{min})^2}{[(T_{opt} - T_{min})(T - T_{opt}) - (T_{opt} - T_{max})(T_{opt} + T_{min} - 2T)](T_{opt} - T_{min})} \end{array} \right.$ |
| No lag Logistic - Cardinal parameter model ⁴ | $\left\{ \begin{array}{l} \frac{dY(t)}{dt} = \mu_{max} [1 - \exp(Y(t) - Y_{max})] \\ \mu_{max} = \frac{\mu_{opt}(T - T_{max})(T - T_{min})^2}{[(T_{opt} - T_{min})(T - T_{opt}) - (T_{opt} - T_{max})(T_{opt} + T_{min} - 2T)](T_{opt} - T_{min})} \end{array} \right.$ |
| Non-isothermal inactivation model | |
| Dynamic Bigelow model ⁵ | $\frac{dy}{dt} = -\frac{1}{D_{ref}} 10^{\frac{T - T_{ref}}{z}}$ |

^{1/2/3} Huang, 2017; ⁵ Van Impe et al., 1992.

9 Statistical indicators in Microrisk Lab

To evaluate and compare the goodness of fit, the statistical indicator of residual sum of squares (RSS, Eq.1), mean square error (MSE, Eq.2), root mean square error (RMSE, Eq.3), regular Akaike information criterion (AIC, Eq.4, Akaike, 1974), modified AIC (AICc, Eq.5, van Boekel et al., 2007; Huang, 2014) and Bayesian information criteria (BIC, Eq.6, Schwarz, 1978). are provided such in the ‘Model Evaluation’ tab.

$$RSS = \sum_{i=1}^n (y_i - \hat{y}_i)^2 \quad \text{Eq.1}$$

$$MSE = \frac{RSS}{n-k} \quad \text{Eq.2}$$

$$RMSE = \sqrt{MSE} \quad \text{Eq.3}$$

$$AIC = n \ln \left(\frac{RSS}{n} \right) + 2(k+1) \quad \text{Eq.4}$$

$$AIC_c = n \ln \left(\frac{RSS}{n} \right) + 2(k+1) + \frac{2(k+1)(k+2)}{n-k-2} \quad \text{Eq.5}$$

$$BIC = n \ln \left(\frac{RSS}{n} \right) + k \ln(n) \quad \text{Eq.6}$$

Where y_i is the i th value of the observation; \hat{y}_i is the i th value of the prediction; k is the number of parameters; n is the number of sample data; L is the maximum value of the likelihood function for the model.

Reference

- Akaike, H., 1974. A new look at the statistical model identification. *IEEE T. Automat. Contr.* 19, 716–723.
- Albert, I., Mafart, P., 2005. A modified Weibull model for bacterial inactivation. *Int. J. Food Microbiol.* 100, 197–211.
- Baranyi, J., Roberts, T.A., 1995. Mathematics of predictive food microbiology. *Int. J. Food Microbiol.* 26, 199–218.
- Bigelow, W.D., 1921. The logarithmic nature of thermal death time curves. *J. Infect. Dis.* 29, 528–536.
- Buchanan, R.L., Golden, M.H., 1995. Model for the non-thermal inactivation of *Listeria monocytogenes* in a reduced oxygen environment. *Food Microbiol.* 12, 203–212.
- Buchanan, R.L., Whiting, R.C., Damert, W.C., 1997. When is simple good enough: a comparison of the Gompertz, Baranyi, and three-phase linear models for fitting bacterial growth curves. *Food Microbiol.* 14, 313–326.
- Geeraerd, A.H., Herremans, C.H., Van Impe, J.F., 2000. Structural model requirements to describe microbial inactivation during a mild heat treatment. *Int. J. Food Microbiol.* 59, 185–209.
- Huang, L., 2014. IPMP 2013—a comprehensive data analysis tool for predictive microbiology. *Int. J. Food Microbiol.* 171, 100–107.
- Huang, L., 2017. Dynamic identification of growth and survival kinetic parameters of microorganisms in foods. *Curr. Opin. Food Sci.* 14, 85–92.
- Huang, L., 2008. Growth Kinetics of *Listeria monocytogenes* in Broth and Beef Frankfurters—Determination of Lag Phase Duration and Exponential Growth Rate under Isothermal Conditions. *J. Food Sci.* 73, E235–E242.
- Huang, L., Hwang, C.-A., Phillips, J., 2011. Evaluating the Effect of Temperature on Microbial Growth Rate—The Ratkowsky and a Bělehrádek-Type Models. *J. Food Sci.* 76, M547–M557.
- Hwang, C.-A., Huang, L., 2018. Dynamic analysis of competitive growth of *Escherichia coli* O157:H7 in raw ground beef. *Food Control* 93, 251–259.
- Koutsoumanis, K.P., Lianou, A., 2013. Stochasticity in Colonial Growth Dynamics of Individual Bacterial Cells. *Appl. Environ. Microbiol.* 79, 2294–2301.
- Presser, K.A., Ratkowsky, D.A., Ross, T., 1997. Modelling the growth rate of *Escherichia coli* as a function of pH and lactic acid concentration. *Appl. Environ. Microbiol.* 63, 2355–2360.
- Ratkowsky, D.A., Lowry, R.K., McMeekin, T.A., Stokes, A.N., Chandler, R.E., 1983. Model for bacterial culture growth rate throughout the entire biokinetic temperature range. *J. Bacteriol.* 154, 1222–1226.
- Rosso, L., Barjard, S., Flandrois, J.P., Lahellec, C., Fournaud, J., Veit, P., 1996. Differential growth of *Listeria monocytogenes* at 4 and 8°C: Consequences for the Shelf Life of Chilled Products. *J. Food Protect.* 59, 944–949.
- Rosso, L., Lobry, J.R., Bajard, S., Flandrois, J.P., 1995. Convenient Model To Describe the Combined Effects of Temperature and pH on Microbial Growth. *Appl. Environ. Microbiol.* 61, 610–616.
- Rosso, L., Lobry, J.R., Flandrois, J.P., 1993. An unexpected correlation between cardinal temperatures of microbial growth highlighted by a new model. *J. Theor. Biol.* 162, 447–463.
- Rosso, L., Robinson, T.P., 2001. A cardinal model to describe the effect of water activity on the growth of moulds. *Int. J. Food Microbiol.* 63, 265–273.
- Schwarz, G., 1978. Estimating the Dimension of a Model. *The Annals of Statistics* 6, 461–464.
- van Boekel, M.A.J.S., 2002. On the use of the Weibull model to describe thermal inactivation of microbial vegetative cells. *Int. J. Food Microbiol.* 74, 139–159.
- van Boekel, M. A. J. S., Zwietering, M. H., 2007. *Experimental design, data processing and model fitting in predictive microbiology*. In: Brul, S., van Gerwen, S., Zwietering, M. (Eds.). Modeling microorganisms in food. Woodhead Publishing, pp. 38.
- Van Impe, J.F., Nicolai, B.M., Martens, T., De Baerdemaeker, J., Vandewalle, J., 1992. Dynamic mathematical model to predict microbial growth and inactivation during food processing. *Appl. Environ. Microbiol.* 58, 2901–2909.
- Vimont, A., Vernozy-Rozand, C., Montet, M.P., Lazizzera, C., Bavai, C., Delignette-Muller, M.L., 2006. Modeling and Predicting the Simultaneous Growth of *Escherichia coli* O157:H7 and Ground Beef Background Microflora for Various Enrichment Protocols. *Appl. Environ. Microbiol.* 72, 261–268.
- Zwietering, M.H., Jongenburger, I., Rombouts, F.M., van 't Riet, K., 1990. Modeling of the bacterial growth curve. *Appl. Environ. Microbiol.* 56, 1875–1881.